

KEN BWD 736 & BWD 738

Large capacity machines for hospitals



Multi purpose

Complete disinfection

Easy handling of utensils

For more information consult - www.ken.dk



Unique advantages with large capacity in minimum space

KEN BWD 736 is available as wall or floor fitted, while BWD738 is a pass-through machine. Even with the door open the BWD is only 88cm in depth. Despite this the machine will hold 3 round bedpans with lid and 3 urine bottles, or various utensils. The bedpan washer can be delivered with manual or automatic door operation.

By automatic door operation it is possible to select between a foot or remote solution. The remote control is placed on a wall so it can be activated by a slight pressure e.g. with the elbow by which means the loading and emptying are more efficient. A wide variety of extra equipment ensures flexibility and optimum use.





Advantages:

- Clear text screen displays information or service messages
- User-friendly operation panel with Proces Indication Bar
- Available with up to 5 programs as standard feature
- Large wash chamber
- Possibility for an utensil basket 430x340x70 mm
- As a standard 1 bedpan and 3 urine bottles
- Pass-through model for isolation rooms

Service

All KEN machines conform to present EU Directives, and all hospital machines have been constructed in accordance with EN ISO 15883, including, for example, the testing of its cleaning capacity according to SIS-TR 3:2001 (formerly SPRI).

Contact your local KEN distributor to make sure you have the correct service for your machines.

The right service will prolong the machine's working life.

Contact a KEN specialist

KEN will gladly help you find the right solution.

Tlf.: +45 62 63 10 91
 e-mail: ken@ken.dk
 Internet: www.ken.dk

KEN partner Australia:
 Malone Hospital Services Pty Ltd.
 24, Tathra Street
 Gosford West 2250, NSW
 Phone: +61 (02) 4323 9731
 E-mail: m-h-s@bigpond.com

Accessories


Look for equipment with the aid of the equipment selector at www.ken.dk, or contact the KEN partner for information about further possible equipment, including customized special solutions.

Multi purpose:

The BWD 736 / 738 can be used for a wide range of ward utensils, e.g. washing bowls, kidney dishes, spitting cups and so on.

Specifications

Power connection (standard)	3N~400/415V, E+Z, 50Hz
Connected load/initial fuse	6kW / 10A alternatively 9kW / 16A
Water connection	Cold water: 1 pcs. ¾" RG, 5-30°C, 0-30° dH
Water connection 2	Hot soft water: 1 pcs. ¾" RG, 30-70°C, 0-3° dH
Water pressure	200 - 800 kPa
Outlet	Ø 110 mm, P/S, wall or floor (to be specified)
Water consumption	
Urine bottle program	2 wash phases, 16 L
Bedpan program	3 wash phases, 24 L
Intensive program	5 wash phases, 40 L
Detergent pump	Standard
Built-in detergent compartment	Room for 2 x 5 L containers
Pump capacity	130 l/min,
Dimension wash chamber HxWxD	500 x 485 x 370mm
Noise level	Approx. 60 dB(A) (sound pressure)
Sealing class	Standard : IP 21 (not scavenging-proof)

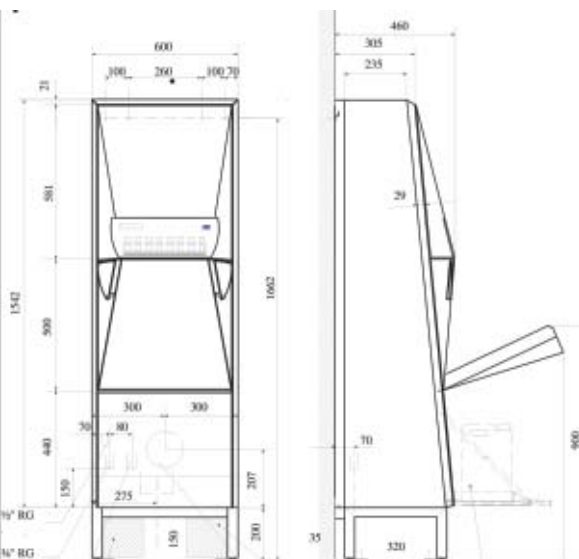

KEN A/S • Art. Nr. 100940 - Technical information and misprint are subject to alterations without notice. 09-2010



We make a difference - locally as well as globally

